

FULLY EQUIPPED FOR ADVENTURE

For people who love to push their boundaries. How much further will you venture?



OFF-ROAD FEATURES & ACCESSORIES

The XT4 offers a range of powered features, seating and accessories designed to support and equip you for adventure. Stay comfortable all day while you enjoy your hobbies, explore new places or hit the beach. With accessories ranging from puncture-resistance kits and camera mounts to sunshades and fishing rod holders, we've thought of all you might need to discover more outdoors.



REHAB OR MPS SEATING

Choose from our rehab or MPS seating and backrests, and a wide range of other options including fold-forward or removable backrests for ease of transport, swing away joystick mounts and a wide range of headrests. All core seating options are standard across the Magic Mobility portfolio.



CHOOSE YOUR BASE COLOUR


Select from our three most popular colours. Our sleek, high-quality metallic finish adds wow factor, and gives depth and dimension to the XT range.



XT4 SPECIFICATIONS

USER WEIGHT LIMIT	with power lift/tilt	160kg / 350lb
	with fixed frame	182kg / 401lb
MAXIMUM SPEED		10km/h / 6.2mph
MAXIMUM SPEED (with lift)		5km/h / 3.1mph
LENGTH OF BASE		1000mm / 39"
WIDTH OF BASE		700mm / 27.5"
WHEELS		Off-road pneumatic 145 / 70-6
SEAT TO FLOOR HEIGHT		480mm / 19"
GROUND CLEARANCE		125mm / 5"
TURNING RADIUS		1410mm / 55.5"
MOTORS		4 x 700W gear in line
CHAIR WEIGHT (minimum)	Base, batteries, rehab seating, centre post legrest, no power lift functions	185kg / 408lb
(standard)	Base, batteries, rehab seating with power lift, tilt, recline and legrest	220kg / 485lb
DRIVE ELECTRONICS		R-net
BATTERY SIZE		115Ah (C20)
SEATING (width & depth)		305-610mm / 12-24"

STANDARD POWER FUNCTIONS

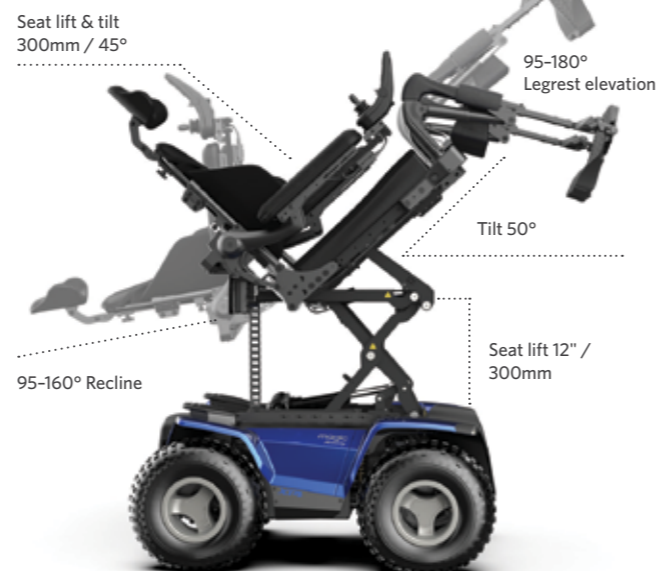

XT4
OFF-ROAD,
OFF THE SCALE



Magic Mobility +61 3 8791 5600
enquiries@magicmobility.com.au
Explore with us magicmobility.com.au



Scan QR to learn more

Please remember that your wheelchair has electric motors and must NEVER be driven through water, into rivers, creeks and the sea. Check online for current version. The information and specifications contained in this brochure are correct and up-to-date at the time of printing and may be subject to change. ©Magic Mobility. Product of Australia.



XT4™

XT4_ENG_SEPT2024_V03

XT4™

VENTURE *Further*

The XT4 is the world's most extreme off-road powerchair in an urban-sized footprint. With front and rear suspension and servo powered steering, it's designed to extend the limits of personal freedom in a world that's not flat.

Delivering more range with high-capacity batteries and integrated high efficiency headlights as standard, the XT4 offers everything you need to venture further - with more confidence - than ever before.



SUPERIOR PATENTED SUSPENSION

Ultimate 4WD traction and stability

Our patented suspension coupled with the power of a motor to each of the four drive wheels delivers the ultimate in obstacle climbing, traction and control.

Front and rear wheels are mounted on their own pivoting arms that adjust independently to follow the contours of the terrain. This helps keep all four wheels on the ground and delivers a smoother ride, greater stability and more control across even the most extreme off-road terrains.

IMPROVED BATTERY RANGE

Go further for longer and see more

Large capacity* batteries come standard on the XT4. These provide exceptional range for a powerful 4WD, extending your scope for adventurous exploration and boosting your confidence to go the distance.

Now you can go further for longer and with added peace of mind.

*115Ah at C20 discharge rate.



PRECISE FORWARD & REVERSE CONTROL

Even on slippery or low traction surfaces, the chair's intuitive steering delivers a feeling of reliability, responsiveness and control, both in forward and reverse. The XT4 steers using a unique combination of differential wheel speeds paired with an industry-leading servo steering module.

INTEGRATED LED LIGHTING

Don't get caught out in the dark! Light your way with confidence thanks to our integrated front and rear LED lighting system.



EATS KERBS FOR BREAKFAST

Experience unstoppable power

With four powerful 700-Watt motors, the XT4 delivers the power you need to plough through practically any terrain and tackle the toughest of life's obstacles. So, if your question is 'Can the XT4 take me there?', the most likely answer is YES!

Because the world's not flat™